

Rabbit Anti MMAE Antibody, mAb

Cat.No: DTA0301


DESCRIPTION

Name	Rabbit Anti MMAE Antibody, mAb
Immunogen	MMAE conjugated to keyhole limpet haemocyanin(KLH)
Host Species	Rabbit
Specificity	This antibody recognizes MMAE and MMAE structurally similar molecules.
Form	Liquid
Concentration	0.5 mg/ml
Storage Buffer	1×PBS, pH7.4, containing 0.1% BSA, 50% glycerol and 0.03% Proclin 300.
Purification	Purified by Protein A
Clonality	Monoclonal
Isotype	Rabbit IgG
Shipping and Storage	Shipped in ice bag. The antibody can be stored for 1-2 weeks at 2-8°C. For long term storage, store it at -20°C. Avoid repeated freezing and thawing cycles.
Expiration Date	Specified on product label.

Applications

Application	Recommended Usage
ELISA	0.1-2 µg/ml
Other applications	User-optimized

Examples

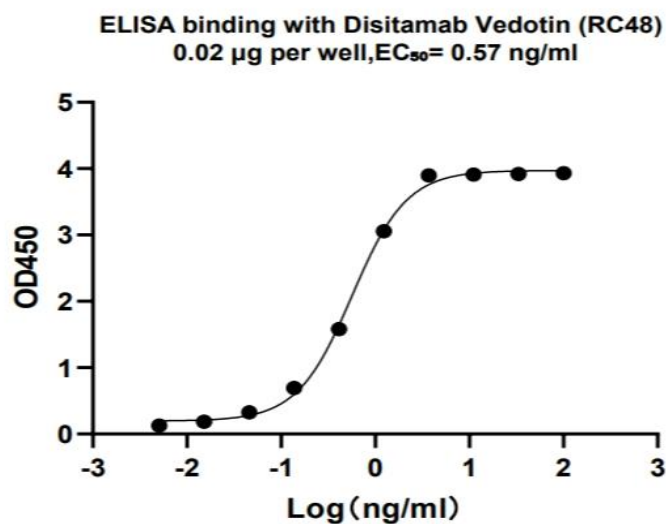

Rabbit Anti MMAE Antibody binds with Disitamab Vedotin (RC48). Coating antigen: Disitamab Vedotin (RC48), 0.2 $\mu\text{g/ml}$, 100 $\mu\text{L/well}$. Rabbit Anti MMAE Antibody dilution start from 100 ng/ml. EC_{50} = 0.57 ng/ml.

感谢您选择德泰生物，我们竭诚为您服务！

For Research Use Only!